



Grib release
master controllers




Confirm
instruments free
1+2




If 3rd arm is free
confirm instrument
free




Surgical gown and
gloves - proceed
directly to OR table

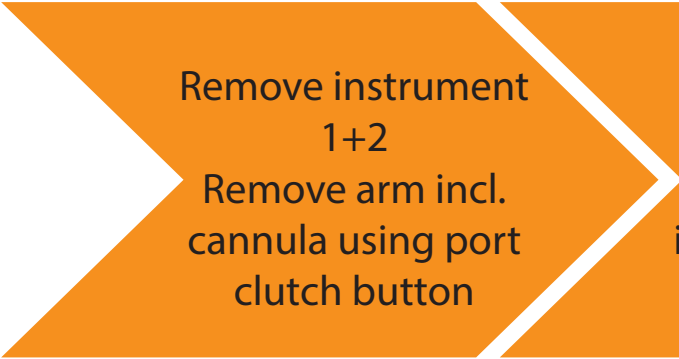

Remove instrument
1+2
Remove arm incl.
cannula using port
clutch button




When instrument 3
is confirmed
released, remove
instrument and arm
incl. cannula

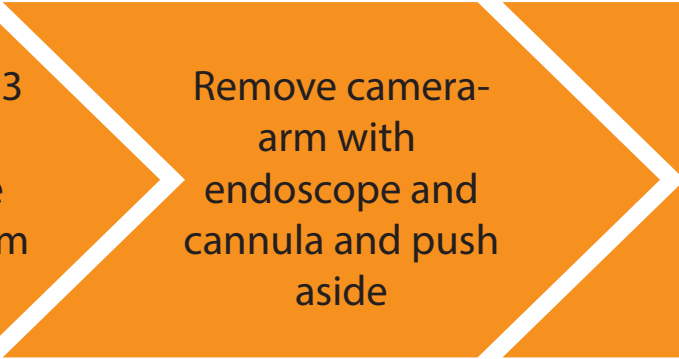

Remove camera-
arm with
endoscope and
cannula and push
aside




NB! If blunt tip is
used, dessufflate




Call for assistance
surgical gown +
groves for console
surgeon



Give instruments
for open procedure
surgeon



Help back up
surgeon with
surgical gown and
gloves

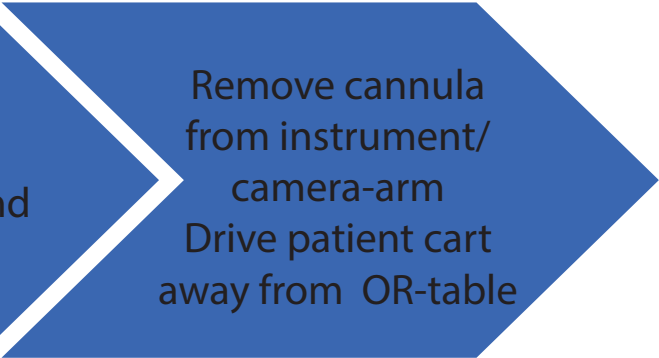

Remove cannula
from instrument/
camera-arm
Drive patient cart
away from OR-table